



Acid Invertase Activity Colorimetric Microplate Assay Kit User Manual

Catalog # CAK1026

(Version 1.3D)

Detection and Quantification of Acid Invertase (AI) Activity in Tissue extracts, Cell lysate Samples.

For research use only. Not for diagnostic or therapeutic procedures.

I. INTRODUCTION.....	2
II. KIT COMPONENTS.....	3
III. MATERIALS REQUIRED BUT NOT PROVIDED.....	4
IV. REAGENT PREPARATION.....	5
V. SAMPLE PREPARATION.....	6
VI. ASSAY PROCEDURE.....	7
VII. CALCULATION.....	8
VIII.TYPICAL DATA.....	9

I. INTRODUCTION

Invertase, also known as sucrase or β -fructofuranosidase, catalyzes the hydrolysis of sucrose by cleaving its glycosidic bond and forming one molecule each of glucose and fructose. According to the optimum pH, invertase is divided into acid invertase (AI) and neutral invertase (NI).

AI mainly exists in cell vacuole or free space, the optimum pH is 4.5-5.0. AI regulates the use of sugar in vacuole and the accumulation of sugar in fruit through the degradation of sucrose in the vacuole.

Acid Invertase Activity Colorimetric Microplate Assay Kit provides a simple and direct procedure for measuring Acid Invertase activity in a variety of samples. AI hydrolyzes the sucrose to generate reducing sugar. The reducing sugar reduces the 3,5-dinitrosalicylic acid to generate red-brown substance. The color intensity, measured at 540 nm, is proportionate to the enzyme activity in the sample.

II.KIT COMPONENTS

Component	Volume	Storage
96-Well Microplate	1 plate	
Assay Buffer	30 mlx 4	4 °C
Reaction Buffer	20 ml x 1	4 °C
Substrate	Powderx 1	4 °C
Dye Reagent	10 mlx 1	4 °C
Standard	Powderx 1	4 °C
Positive Control	Powderx 1	-20 °C
Plate Adhesive Strips	3 Strips	
Technical Manual	1 Manual	

III. MATERIALS REQUIRED BUT NOT PROVIDED

1. Microplate reader to read absorbance at 540 nm
2. Distilled water
3. Pipettor, multi-channel pipettor
4. Pipette tips
5. Mortar
6. Ice
7. Centrifuge
8. Timer
9. Convection oven

IV. REAGENT PREPARATION

Standard: Briefly centrifuge prior to opening. Dissolve in 1 ml distilled water to generate 10 mmol/L of standard stock solution. Store at 4 °C for 1 month or -20°C for 6-7 months after reconstitution. Then dilute to 3 mmol/L standard top solution by adding 0.3 ml stock solution into 0.7 ml distilled water. Perform 2-fold serial dilutions of the top standard solution using distilled water to make the standard curve. The concentration of standard curve could be 3/1.5/0.75/0.375/0.188/0.094/0.047 mmol/L.

Substrate: Add 8 ml Reaction Buffer to dissolve before use. Store at 4 °C for 3-4 months after reconstitution.

Positive Control: Briefly centrifuge prior to opening. Add 1 ml Assay Buffer to dissolve before use, then add 0.1 ml into Assay Buffer 0.9 ml. Store at -80°C for 1 month after reconstitution.

Note: Divide into small aliquots to avoid repeated freeze-thaw cycles.

V. SAMPLE PREPARATION

1. For tissue samples

Weigh out 0.1 g tissue, homogenize with 1 ml Assay Buffer on ice, centrifuged at 12000g 4°C for 10 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

2. For cell and bacteria samples

Collect cell or bacteria into centrifuge tube, discard the supernatant after centrifugation, add 1 ml Assay Buffer for 5×10^6 cell or bacteria, sonicate (with power 20%, sonicate 3s, interval 10s, repeat 30 times); centrifuged at 12000g 4°C for 10 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

VI. ASSAY PROCEDURE

Add following reagents in the microcentrifuge tubes:

Reagent*	Sample**	Control	Standard	Blank	Positive Control
Sample	20 μ l	--	--	--	--
Reaction Buffer	--	20 μ l	--	--	--
Positive Control	--	--	--	--	20 μ l
Substrate	80 μ l	80 μ l	--	--	80 μ l
Mix, put it in the oven, 37°C for 10 minutes.					
Standard	--	--	100 μ l	--	--
Distilled water	--	--	--	100 μ l	--
Dye Reagent	100 μ l	100 μ l	100 μ l	100 μ l	100 μ l
Mix, put it into the convection oven, 90 °C for 10 minutes, record absorbance measured at 540nm.					

Note:

*Reagents must be added sequentially and should not be premixed prior to addition.

**For unknown samples, we recommend doing a pilot experiment & testing several doses to ensure the readings are within the standard curve range. If the enzyme activity is lower, please add more samples into the reaction system; or increase the reaction time; if the enzyme activity is higher, please dilute the sample, or decrease the reaction time.

**For colored samples, we recommend setting a parallel sample background control well with same volume of sample only. Other reagents were replaced by distilled water to the same total volume. Subtract the OD value of the sample background control from the OD value of the sample to correct for interference from the sample's own color.

VII. CALCULATION

Unit Definition: One unit of AI activity is defined as the enzyme generates 1 μmol of reducing sugar per minute.

1. According to the slope of the standard curve

$$\text{Activity} = \frac{(\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) - \text{Intercept}}{\text{Slope} \times T} \times \frac{V_{\text{Standard}}}{V_{\text{Sample}}} \times n (\text{U/mL})$$

2. According to one point of the standard OD value and concentration

2.1 According to the protein concentration of sample

$$\text{Activity} = \frac{(C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}})}{(\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \times V_{\text{Sample}} \times C_{\text{Protein}} \times T} (\text{U/mg})$$

2.2 According to the quantity of cells or bacteria

$$\text{Activity} = \frac{(C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}})}{(\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \times N \times (V_{\text{Sample}} / V_{\text{Assay}}) \times T} (\text{U}/10^4)$$

2.3 According to the weight of sample

$$\text{Activity} = \frac{(C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}})}{(\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \times W \times (V_{\text{Sample}} / V_{\text{Assay}}) \times T} (\text{U/g})$$

Slope: the absorbance slope of standard curve

n: the dilution factor

C_{Protein} : the protein concentration of sample, mg/mL

W: the weight of total sample, g

N: the quantity of total cell or bacteria sample, 10^4

C_{Standard} : the concentration of standard, mmol/L = $\mu\text{mol/mL}$

V_{Standard} : the volume of standard in assay procedure, mL


V_{Sample} : the volume of sample in assay procedure, mL

V_{Assay} : the volume of Assay Buffer in sample preparation, mL


T: the reaction time, minute

VIII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

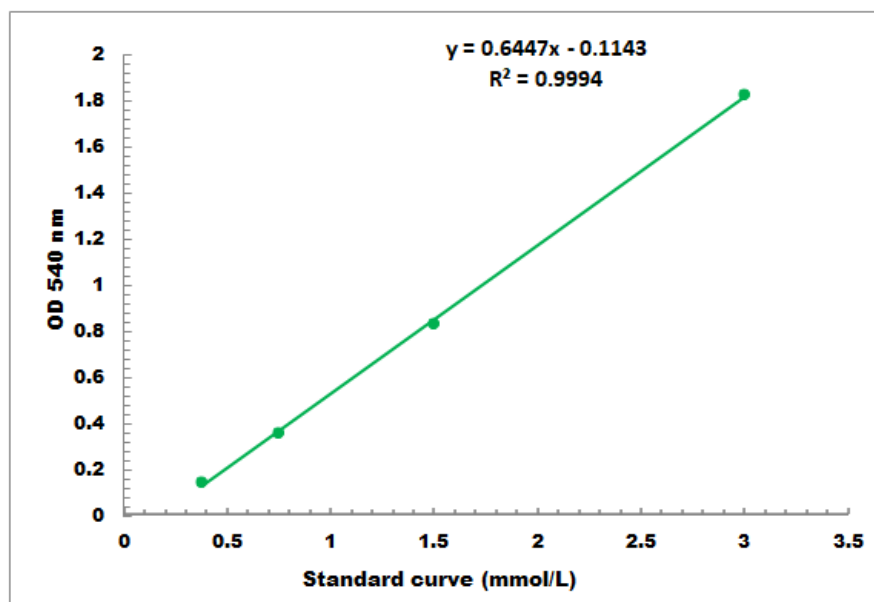

Detection Range: 0.3mmol/L -3mmol/L

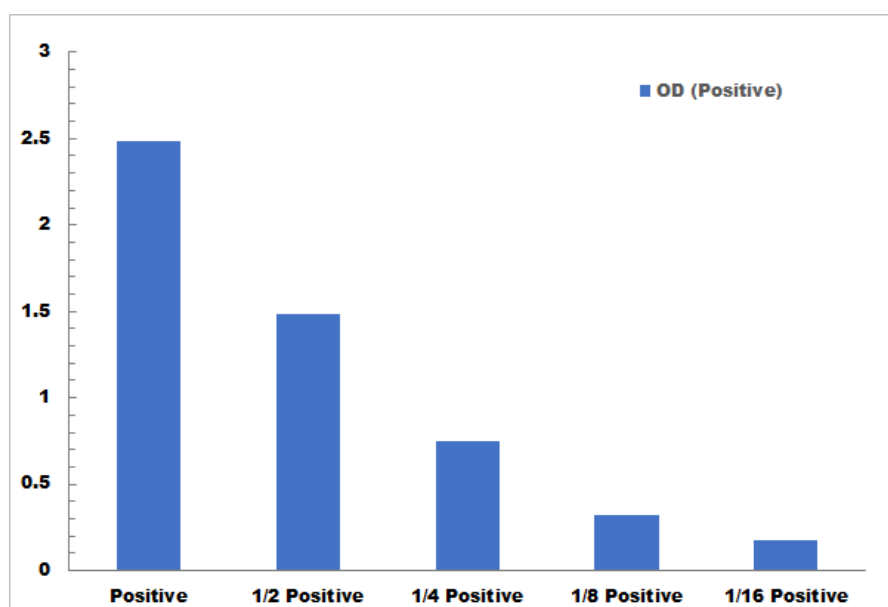

Positive Control reaction in 96-well plate assay with decreasing the concentration