


Product Data Sheet

Human ATP6V1H Knockdown HEK293T Cell Lysate


Catalog #	Source	Reactivity	Applications
CKD0450	HEK293T	Human	
Description	Human ATP6V1H knockdown HEK293T cell lysate is engineered by our optimized transduction of the specific shRNA.		
Form	Liquid in 50 mM Tris-HCl, 2% SDS, 10% glycerol, 100 mM DTT, 0.01% Bromophenol blue		
Gene Symbol	ATP6V1H		
Alternative Names	V-type proton ATPase subunit H; V-ATPase subunit H; Nef-binding protein 1; NBP1; Protein VMA13 homolog; V-ATPase 50/57 kDa subunits; Vacuolar proton pump subunit H; Vacuolar proton pump subunit SFD		
Entrez Gene	51606		
SwissProt	Q9UI12		
KO/KD Validation	RT-qPCR		
Cell Line	HEK293T		
Components	Human ATP6V1H knockdown HEK293T cell lysate x 1 vial; shControl HEK293T cell lysate x 1 vial		
Directions for Use	Thaw the vial on ice or at room temperature. Mix gently by pipetting or brief vortexing. Briefly centrifuge to collect the sample. Load the desired amount of lysate directly onto an SDS-PAGE gel. Perform SDS-PAGE and Western blot analysis using standard laboratory protocols. Note: This lysate has been pre-denatured in SDS sample buffer containing DTT and does not require additional reducing agent or boiling before electrophoresis.		
Storage/Stability	Shipped at 4°C. Upon delivery aliquot and store at -80°C for one year. Avoid freeze/thaw cycles.		

COHESION BIOSCIENCES LIMITED


Product Data Sheet

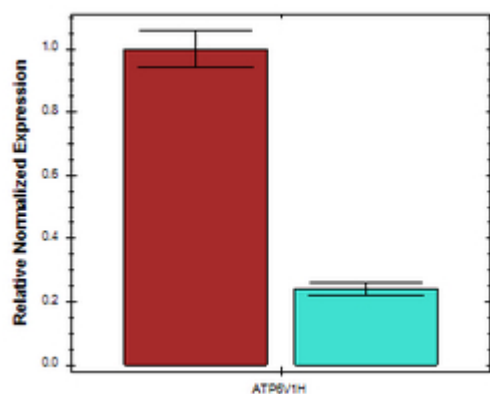


RT-qPCR analysis of ATP6V1H knockdown efficiency. Compared with the 293T shRNA control (left), ATP6V1H mRNA expression in the stable shRNA knockdown cell line (right) decreased, and the knockdown efficiency was 76%.

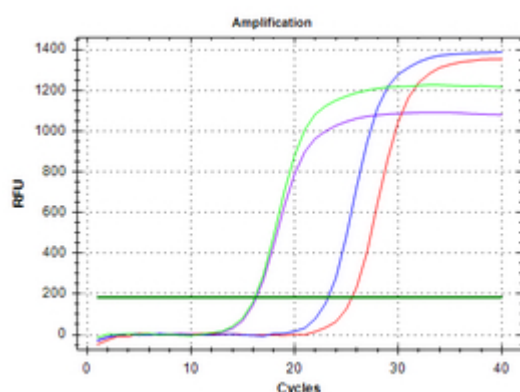

RT-qPCR amplification curves. The Ct values of the knockout sample significantly increased, confirming the effective knockdown of ATP6V1H mRNA. [Red curve (knockdown sample), Blue curve (shRNA control), Purple and Green curve (reference gene GAPDH)]

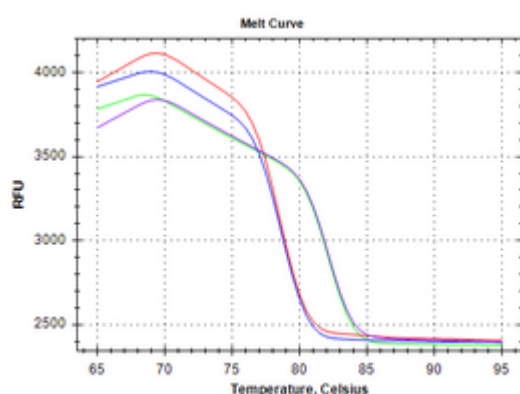

RT-qPCR melt curve analysis demonstrated highly specific PCR amplification product without detectable nonspecific amplification or primer-dimer formation. [Red curve (knockdown sample), Blue curve (shRNA control), Purple and Green curve (reference gene GAPDH)]

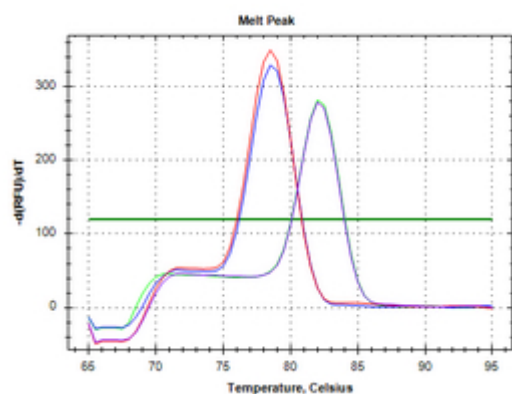

This picture presents the raw fluorescence data (relative fluorescence units, RFU) used to generate the melt curves. [Red curve (knockdown sample), Blue curve (shRNA control), Purple and Green curve (reference gene GAPDH)]

COHESION BIOSCIENCES LIMITED