

Product Data Sheet


Human SGMS2 Knockdown HEK293T Cell Lysate

Catalog #	Source	Reactivity	Applications
CKD0509	HEK293T	Human	
Description	Human SGMS2 knockdown HEK293T cell lysate is engineered by our optimized transduction of the specific shRNA.		
Form	Liquid in 50 mM Tris-HCl, 2% SDS, 10% glycerol, 100 mM DTT, 0.01% Bromophenol blue		
Gene Symbol	SGMS2		
Alternative Names	SMS2; Phosphatidylcholine:ceramide cholinephosphotransferase 2; Sphingomyelin synthase 2		
Entrez Gene	166929		
SwissProt	Q8NHU3		
KO/KD Validation	RT-qPCR		
Cell Line	HEK293T		
Components	Human SGMS2 knockdown HEK293T cell lysate x 1 vial; shControl HEK293T cell lysate x 1 vial		
Directions for Use	Thaw the vial on ice or at room temperature. Mix gently by pipetting or brief vortexing. Briefly centrifuge to collect the sample. Load the desired amount of lysate directly onto an SDS-PAGE gel. Perform SDS-PAGE and Western blot analysis using standard laboratory protocols. Note: This lysate has been pre-denatured in SDS sample buffer containing DTT and does not require additional reducing agent or boiling before electrophoresis.		
Storage/Stability	Shipped at 4°C. Upon delivery aliquot and store at -80°C for one year. Avoid freeze/thaw cycles.		


COHESION BIOSCIENCES LIMITED

WEB
www.cohesionbio.comORDER
order@cohesionbio.comSUPPORT
techsupport@cohesionbio.comCUSTOM
custom@cohesionbio.com


Product Data Sheet

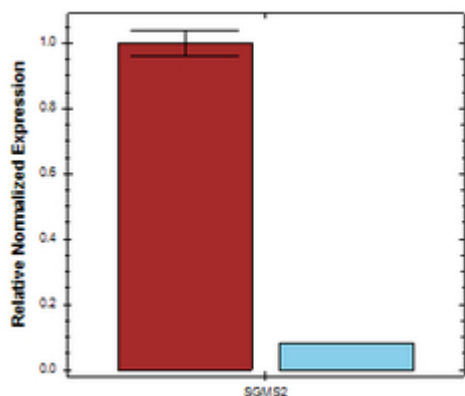


RT-qPCR analysis of SGMS2 knockdown efficiency. Compared with the 293T shRNA control (left), SGMS2 mRNA expression in the stable shRNA knockdown cell line (right) decreased, and the knockdown efficiency was 92%.

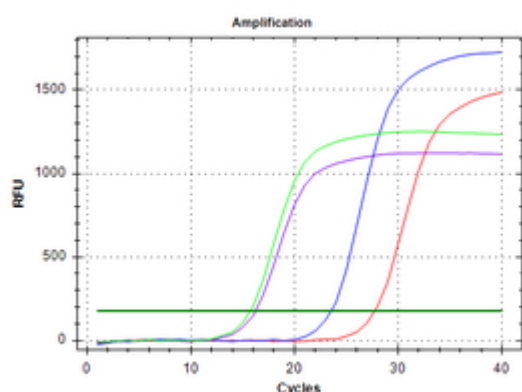

RT-qPCR amplification curves. The Ct values of the knockout sample significantly increased, confirming the effective knockdown of SGMS2 mRNA. [Red curve (knockdown sample), Blue curve (shRNA control), Purple and Green curve (reference gene GAPDH)]

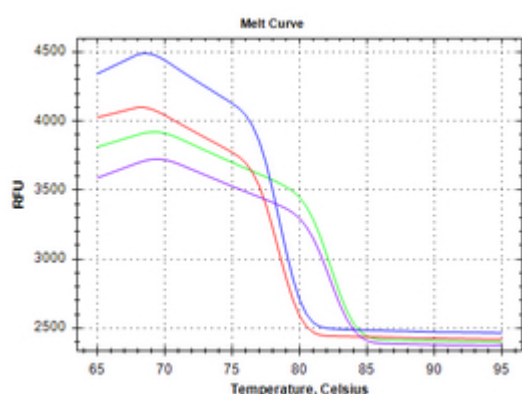

RT-qPCR melt curve analysis demonstrated highly specific PCR amplification product without detectable nonspecific amplification or primer-dimer formation. [Red curve (knockdown sample), Blue curve (shRNA control), Purple and Green curve (reference gene GAPDH)]

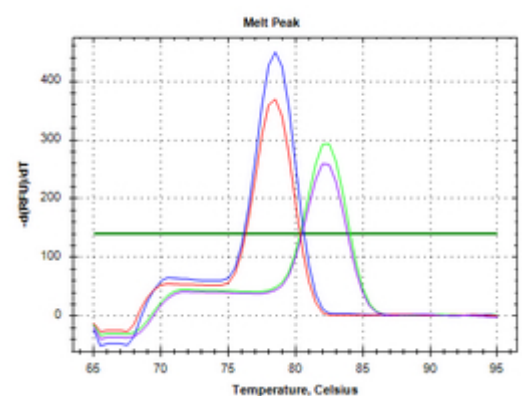

This picture presents the raw fluorescence data (relative fluorescence units, RFU) used to generate the melt curves. [Red curve (knockdown sample), Blue curve (shRNA control), Purple and Green curve (reference gene GAPDH)]

COHESION BIOSCIENCES LIMITED